

**1**

**Blood Culture Bottles**  
 (Blood Culture Kit or Pediatric)

*\*If using a butterfly needle, be sure to draw a spare tube before drawing the Sodium Citrate (Light Blue Top) Tube. This can be an Serum Separator Tube (SST), Tiger Top, or Red Top Tube (No Gel)*



**2**

**Sodium Citrate (Light Blue Top) Tube**

*\*Must be filled to fill line on tube*



**3**

**Red Top Tube (No Gel)**



**4**

**Serum Separator Tube (SST)**

**Tiger Top Tube**

**Royal Blue Trace Element Serum Tube**



**5**

**Lithium Heparin Tube**

**Sodium Heparin Tube**

**Royal Blue Trace Element Sodium Heparin Tube**





**6**

**EDTA Lavender Top Tube**

**Plasma Preparation Tube (PPT)**

**Royal Blue Trace Element K2 EDTA Tube**

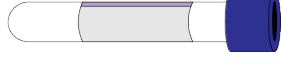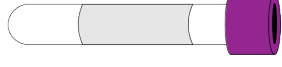

**7**

**Sodium Fluoride/Potassium Oxalate Gray Top Tube**



**8**

**Yellow ACD A/B Solution Tube**

